


Specifications of the High Standard series of Dhona Balances

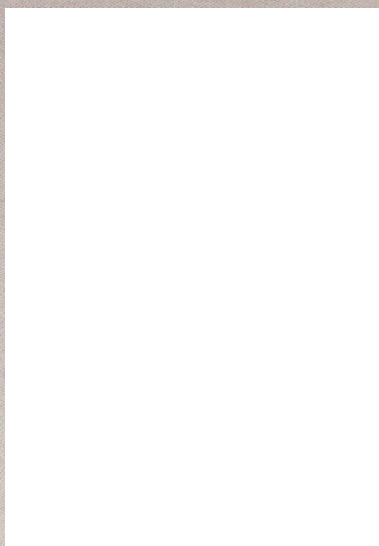

MODELS

Readability

Capacity

Optical Range


Scale Division

Built-in Mechanical weight

Accuracy of Optical Range

Weighing Principle

Power Connection

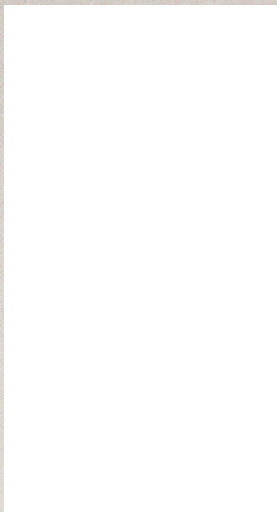

200 D

0.1 mg.

200 gm.

1000 mg.

10 mg.

1-199 gm.

± 0.1 mg.

Substitution

220/230 Volt. AC.

100 DS

0.01 mg.

100 gm.

100 mg.

1 mg.

0.1-99.9 gm.

± 0.01 mg.

Substitution

220/230 Volt. AC

Standard features of "DHONA" Analytical Balance :

- Digital (Mechanical) indication with large figures for rapid and precise reading
- Clearly-lit scale window - the figures immediately above and below the measured value are also visible through green transparent shades.
- Very short settling time and high weighing accuracy.
- All controls conveniently located for efficient and reliable operation
- Drop weight system with built-in drop weights in anti-magnetic stainless steel.
- Zeroing device for setting the loaded machine to zero reading.
- Stainless Steel pan.
- Weighing Space with sliding side doors and fixed glass front.
- Floor of weighing space also of glass.
- Powder coat-painted pressure die cast Aluminium Housing